



UNIVERSAL SCREEN FOR ALL REFRACTION ROOMS



A large LCD screen which adapts to all refraction rooms. Its design and great possibilities of customisation, are highly appreciated.







CS pola 600

CHART SYSTEM







pola





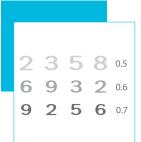



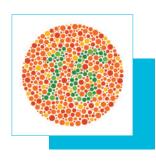



BRIGHT & CLEAR IMAGES

• 24" screen offering wide viewing angle to provide clear chart and remove image distortion.

WIDE RANGE OF TESTS

- Visual acuity: Landolt C, numbers, letters, including ETDRS.
- Special tests: binocular and fusion, aniseikonia, stereopsis, strabismus, phoria, color blindness, Amsler grid test...
- Prescription check in real life situation.
- Random function to avoid patient optotypes memorisation.

EASY OPERATION

- Remote control with built-in LCD panel for optotypes visualisation.
- Can be also connected to the APH 500 to build a fully automated refraction chain.

| SCREEN SIZE | 24" | TEST CHARTS | |
|---------------------|-------------------------|-----------------------------|--|
| LUMINOSITY | Max 300Cd/M² | Optotypes | Landolt C, Numbers, Letters, ETDRS |
| RESOLUTION | 1920 x 1080 pixels | | Astigmatism, Cross-grid, Cross cylinder Lines (Horizontal/Vertical) Duochrome balance - Binocular balance Fixation point - Aniseikonia Stereopsis - Phoria Ishihara test - Amsler grid test Red/Green - Worth. |
| PROJECTION DISTANCE | 1.5 to 8m (0.2m) | | |
| PROJECTION SETTINGS | - | Masks and tests | |
| CHART WINDOW SIZE | 298.89(H) x 531.36(W)mm | | |
| | | DIMENSIONS AND WEIGHT | 460(W) x 392(H) x 65(D)mm - 7.5kg |
| | | Power supply | AV 100-120V/220-240V, 50-60Hz |
| | | Power consumption | Less than 80W |
| | | C Conformity marking | |

PRODUCT FUNCTIONS











Instruments Division 64 bis, avenue Aubert 94306 Vincennes Cedex France Tel: 33 (0)1 72 94 71 00

Fax: 33 (0)1 72 94 70 27 www.essilor-instruments.com

Your local distributor: